


### Features:

The Axia SB4 is a four-movement-imaging table (up/down, lateral tilt, trendelenburg and longitudinal travel) intended to be used with mobile x-ray units and c-arms for minimally invasive procedures. This imaging table comes with standard 500lbs load rating, 2" patient pad for patient comfort, 24" wide top, locking swivel castors, narrow base for easy patient transfer, simultaneous motion movements and extremely low boarding height. Axia SB4' diving board style is designed for easy c-arm access. The tabletop material is made up of carbon fiber skin/foam core with an integral facial cut-out with filler piece

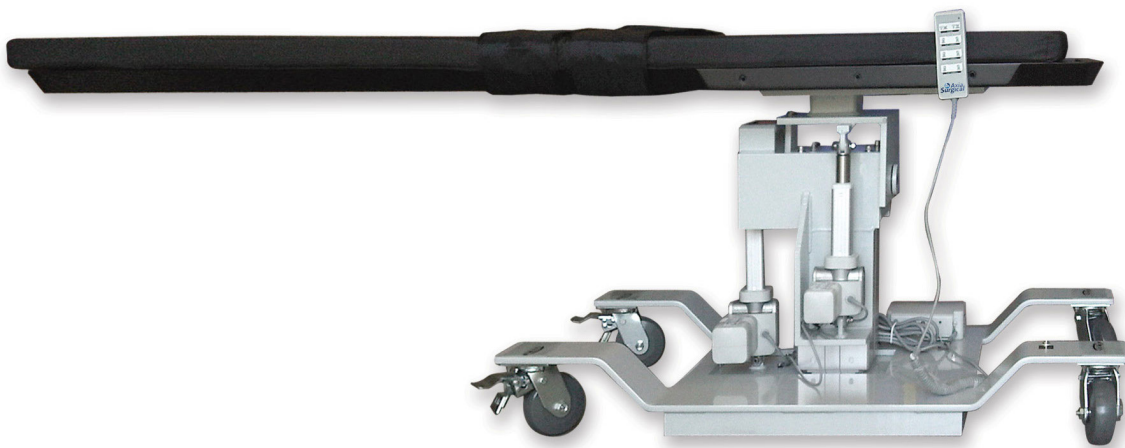

### Specifications

- **Height:** 29" ~ 39"
- **Trendelenburg:** (+/-) 13°
- **Lateral Tilt:** (+/-) 17°
- **Longitudinal travel:** 10" (+/-5")
- **Patient Capacity:** 500lbs (227kg)
- **Tabletop and Mattress Attenuation:** 1.2mm Al. (+/- .3mm tolerance)
- **Radiolucent Area:** 56"L x 24"W unobstructed
- **Tabletop:** 84"L x 24"W x 1.62"T
- **Table Dimensions:** 86"L x 26"W
- **Table Weight:** 740lbs (335kg)